


MRI DISTORTION




KEY FEATURES

- The **GRID** Phantom is designed for MR Image Distortion assessment
- Verify MR Image Distortion in 3 Axes to Identify Hardware-Related Geometric Inaccuracies
- The 3D volume of the **GRID** Phantom is comprised of 9 Image Planes with 4,689 individual Spherical Fiducials
- Maximum Scan Plane Field of View
- Can be scanned in Axial, Saggital, and Coronal orientations
- Solid Acrylic Construction

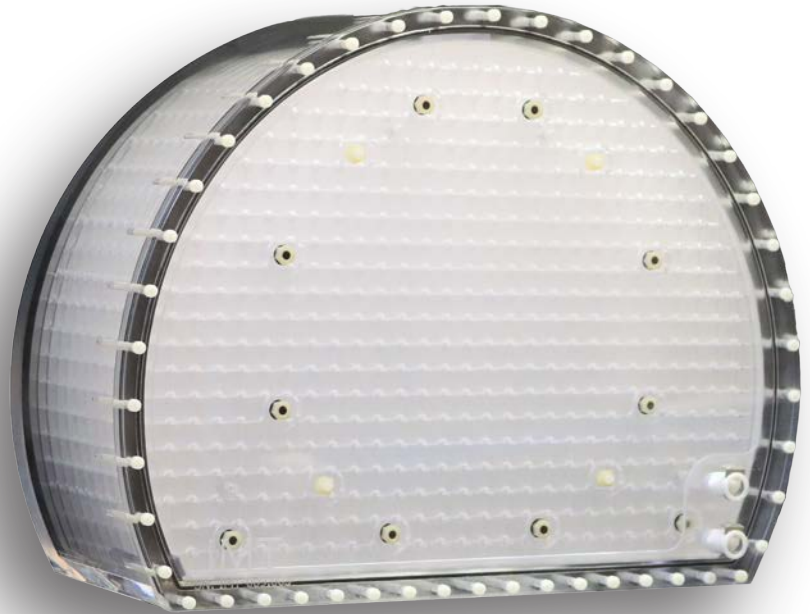


FEATURING **iQA**
AUTO ANALYSIS SOFTWARE

SPATIAL ANALYSIS

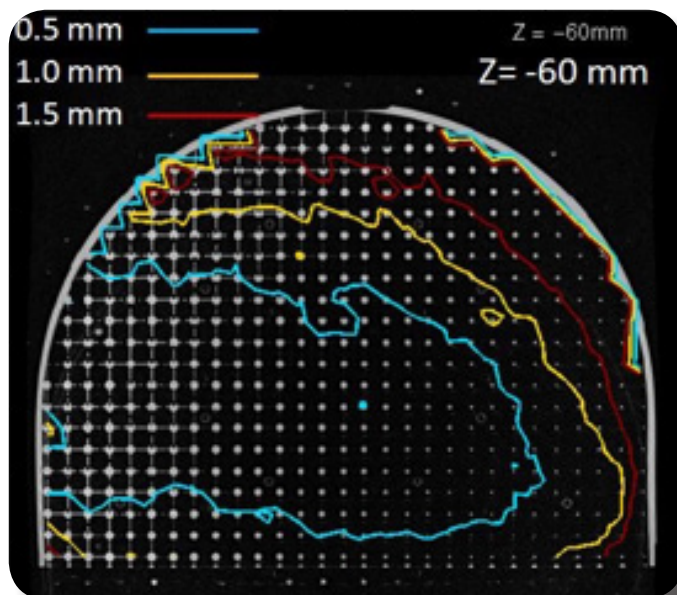

AXIAL IMAGE

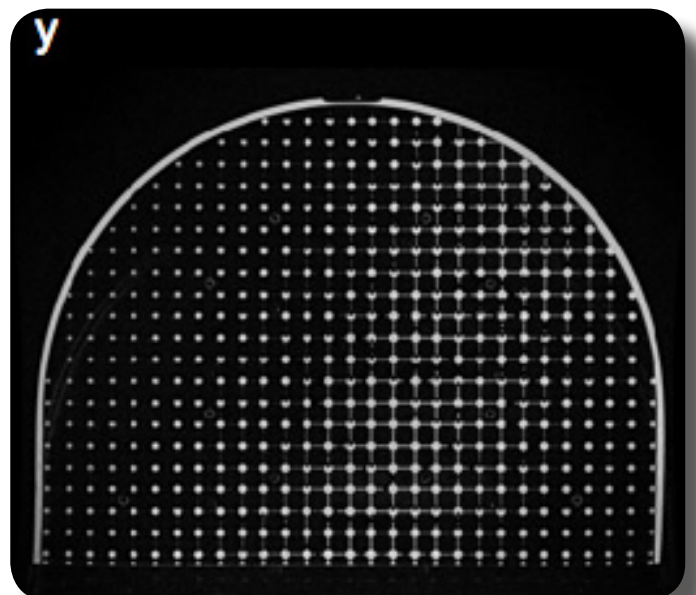

IMT



深圳为尔康科技有限公司

联系人：曾祥满 手机：13632925349 QQ：274798107

电话：0755-28896837 地址：深圳市龙岗区沙平北路111号608A

网址：www.medicalQC.com 邮箱：szchina1718@163.com

MRI DISTORTION



9 IMAGE PLANES



SAGGITAL SCAN ORIENTATION



DISTORTION MEASUREMENT

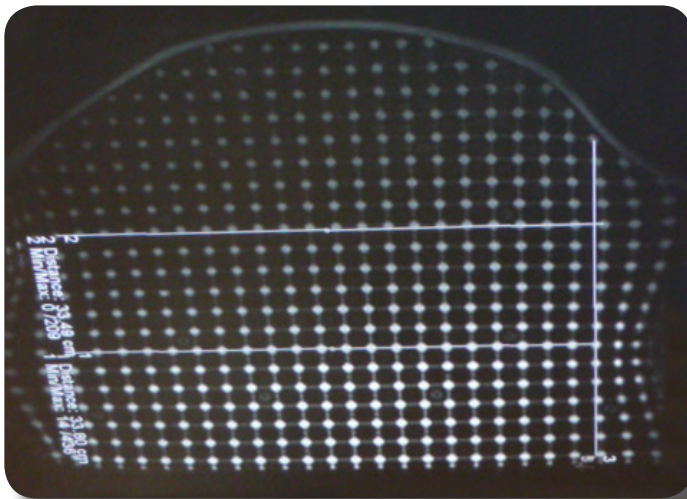


IMAGE DISTORTION




GRID™ SPECIFICATIONS


PHYSICAL SIZE: 465 x 350 x 168 mm

WEIGHT (FILLED): 40 Kg

D = 9mm, Center-to-Center 16mm (+/- 0.1mm)

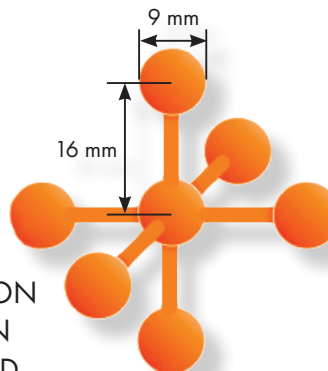

CARRY STRAP INCLUDED

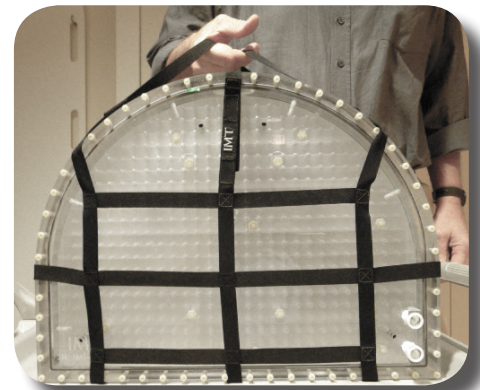

IMT

NOT FOR DISTRIBUTION
OR REPRODUCTION
WITHOUT EXPRESSED
CONSENT