


DR QC Phantom

Model 139701

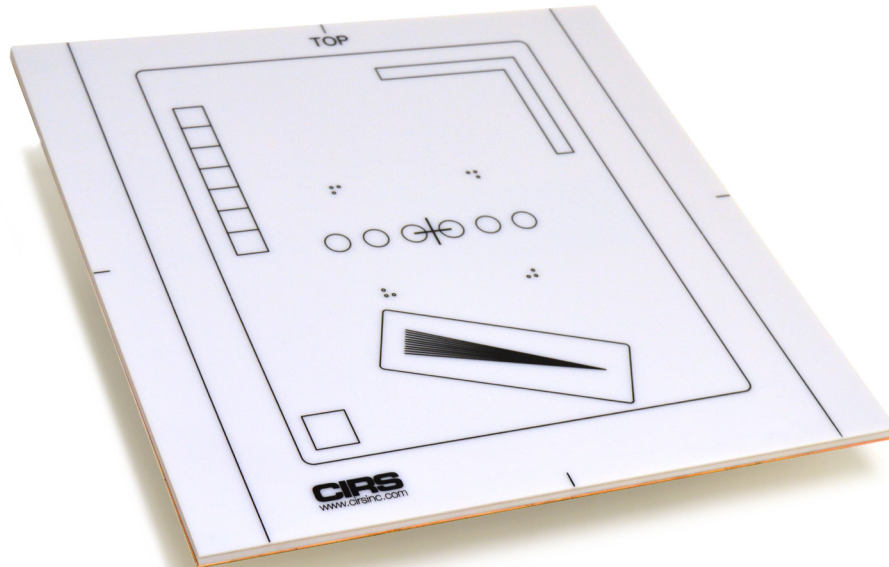

QUALITY ASSURANCE FOR DR SYSTEMS

The CIRS DR QC Phantom allows comprehensive quality assurance testing for DR Detectors and associated software.

The phantom enables seven quality assurance assessments for monthly and semi-annual evaluation, including exposure linearity and sensitivity, high and low contrast reproducibility, artifact and residual image detection, image resolution, collimator beam alignment, measurement tool accuracy, and display jitter QA testing.

The CIRS DR QC Phantom measures 17" x 17" x 0.5" and allows for quick checks on 14" x 17" or 17" x 17" DR systems. Lead content has been minimized to meet European regulations.

The phantom comes with standard suggested testing parameters, QC report forms, soft carry case and 48-month warranty.

Features

- Allows testing of 14"x17" (35 x 43 cm) and 17"x17" (35 x 35 cm) DR systems
- Includes test objects for seven quality assurance tests
- Complies with German standards DIN 6868-58 and DIN 6868-13
- Contains DIN Type 23 Linepair target
- Minimized lead content to meet European regulations
- No moving parts
- Testing parameters and QC Report forms provided
- Extended 48-month warranty


* This product was developed in collaboration with Del Medical of Bloomingdale, Illinois.



深圳为尔康科技有限公司

曾生：13632925349 QQ：274798107 电话：0755-28896837

地址：深圳市龙岗区沙平北路111号吉茂大厦608A 网址：www.medicalqc.com

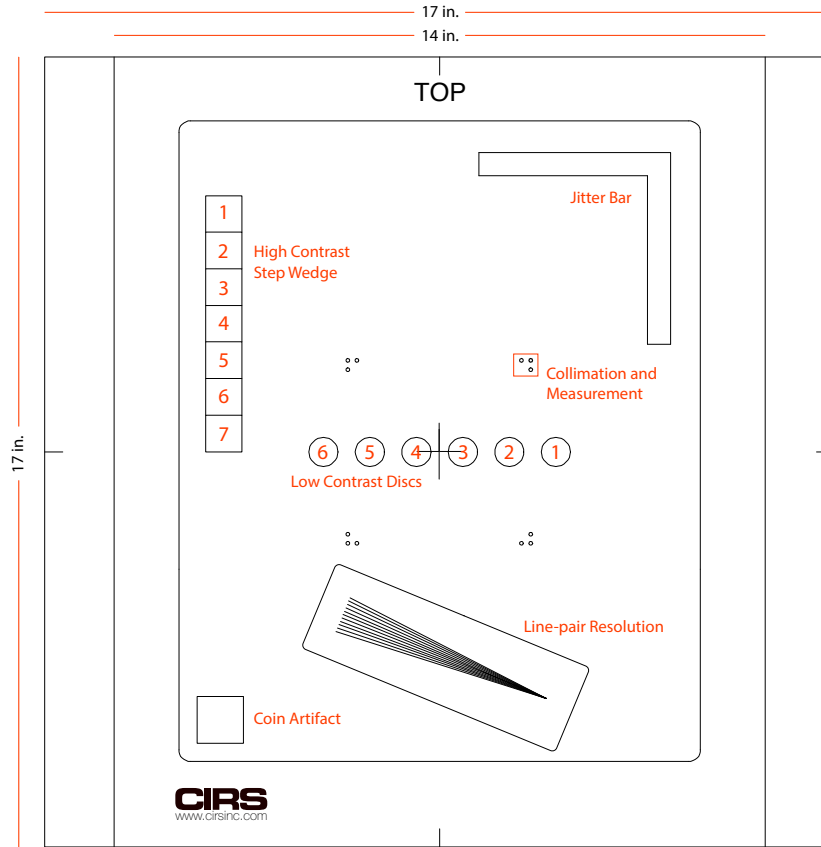

Suggested Test Frequency

Monthly Testing

- Exposure linearity and sensitivity
- Tracking contrast reproducibility high and low contrast
- Checking for artifacts and residual images
- Image Resolution

Semi Annual Testing

- Collimator Beam Alignment
- Accuracy of measurement tools
- Display Jitter evaluation

SPECIFICATIONS

DIMENSIONS:	17"x17"x~1/2"
WEIGHT:	8 lb
MATERIALS:	PMMA, Copper, Aluminum

MODEL 139701 INCLUDES

QTY	COMPONENT DESCRIPTION
1	DR QC Phantom
1	Padded Soft Case
1	User Guide
-	48 Month Warranty