DEXA Phantom

Model 026



THE "BONA FIDE PHANTOM" (BFP)(1)

BioClinica's "Bona Fide Phantom" (BFP), is a quality control tool for Dual-Energy X-ray Absorptiometry (DEXA) instruments, which features an acrylic-embedded calcium hydroxyapatite (CHA) step-wedge. Advanced design features make it the best choice for assessing DEXA instrument stability. You can successfully use the BFP on all mainstream DEXA instruments.

The BFP offers a range of densities (0.7 - 1.5 g/cm2), to verify instrument function over the clinically relevant range, not just at a single "healthy" BMD. Linearity of BMD over this range is critical for full instrument evaluation. Our phantom uses a CHA insert for direct assessment of bone density accuracy, as required for cross-calibration phantoms used in clinical trials.

The BFP is cast in acrylic and comes with its own carrying case for easy handling. The tote remains on the phantom during scanning and does not affect BMD readings, allowing rapid placement and removal for the phantom from the bed. A flight case is available as an option.

Features

- Checks a Range of Densities
- Easy to carry and can be Scanned in the Bag
- Universal Axial DEXA Instrument Compatibility
- No Water Bath
- Edge-Detection Challenge
- Realistic Soft Tissue Mimic


(1) BFP design is the property of BioClinica, Inc.



深圳为尔康科技有限公司 联系人:曾祥满 手机:13632925349
QQ:274798107 电话:0755-28896837 地址:深圳市龙岗区沙平北路111号6008

网址:www.medicalQC.com 邮箱:szchina1718@163.com

DEXA PHANTOM Model 026

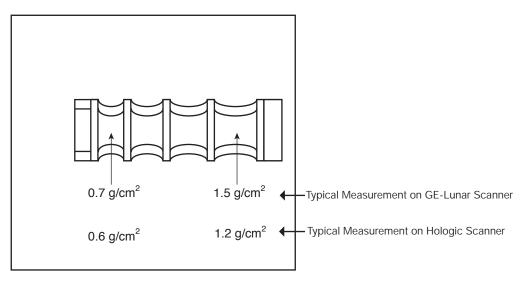

Note: The various DXA scanner manufacturers have developed and published cross-calibration formulas for use in data comparison.

SPECIFICATIONS

OVERALL DIMENCIONS	
OVERALL DIMENSIONS	22 cm x 19 cm x 15 cm
	(8.7" x 7.5" x 6")
WEIGHT	17 lbs. (6.3 kg)
MATERIALS	ACRYLIC EMBEDDED Calcium Hydroxyapatite
	SOFT TISSUE ANALOG Approximately 25 % Fat
	BMD RANGE 0.7 to 1.5 g/cm² (GE-Lunar) 0.6 to 1.2 g/cm² (Hologic)

MODEL 026 INCLUDES

QTY	COMPONENT DESCRIPTION
1	DEXA Phantom
1	Carrying Case
1	User Guide
-	60-Month Warranty



深圳为尔康科技有限公司 联系人:曾祥满 手机:13632925349

QQ: 274798107 电话: 0755-28896837 地址: 深圳市龙岗区沙平北路111号6008

网址:www.medicalQC.com 邮箱:szchina1718@163.com